



People.Health.Care.

Vliwasorb® sensitive border **NEW**

Easy application,
good exudate management,
reduction of germs.



Vliwasorb® sensitive border **NEW**

Superabsorbent wound dressing with silicone and adhesive border.

Polyurethane film layer

- high MVTR
- waterproof & with barrier against external contaminants

Breathable hydrophilic enclosed layer

- surrounding the absorbent core

Superabsorbent core

- for absorbing and retaining exudate, which reduces the risk of maceration
- wound exudate may contain germs, exudate containing germs is bound within the dressing and safely removed during dressing change

Liquid permeable layer

- surrounding the absorbent core
- the exudate is evenly distributed

Silicone wound contact layer

- rapid absorption and transfer of the exudate to the absorbent core
- larger pores enable the absorption of viscous exudates
- gentle to the skin

Indications:


Degree of exudation



For moderately to heavily exuding wounds


Wound type

Acute wounds

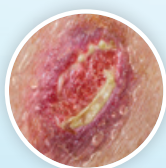


Surgical wounds


Chronic wounds

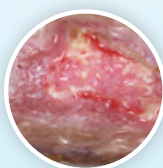

Leg ulcers

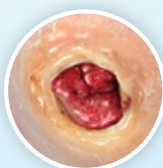

Diabetic foot ulcers

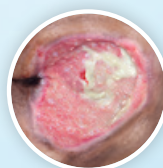

Pressure ulcers

Wound depth



Superficial wounds



Deep wounds*

Vliwasorb sensitive border can ...

- be used under compression.
- be retained with a suitable material (e.g. Curafix H, Mollelast).
- do not cut the wound pad. The border can be cut into a concave shape or cut apart. ✗

* The product can be used as a secondary dressing treating deep wounds when combined with an appropriate wound filler (e.g. Suprasorb A (Pro), Suprasorb A+Ag (Pro), Suprasorb Liquacel Pro, Suprasorb C).

Vliwasorb® sensitive border NEW

Your benefits.

The product features

Your benefit



Easy application and atraumatic removal

- atraumatic and virtually painless removal/dressing changes
- **user-friendly adherent edge**, which facilitates the attachment of the dressing and allows repositioning



Good exudate management

- the high absorption capacity allows long dressing change intervals and long periods of wound rest
- larger pores enable **the absorption of viscous exudates**
- **reduces the risk of maceration**



Reduction of germs

- reduction of germs by **binding in the dressing** (to the superabsorbent polymers); **removal during dressing changes**

Vliwasorb® sensitive border ...

... retains its shape!

Good shape retention even when used in a vertical position as a result of **compressed superabsorbent polymer layer**.

Vliwasorb
sensitive
border
before use



Vliwasorb
sensitive
border
with maximum
absorption

... is practical to use and good for the patient:

- it is gentle to the skin
- offers gentle comfort, minimises pain during removal and reduces trauma to the wound and the wound environment
- clothing and bed linen are protected from soaking by the practical clothing protection layer

... has the ability to absorb and lock in exudate.

Wound exudate may contain **bacteria, fungi** (such as *Candida albicans*, *Candida parapsilosis*) and **endotoxins**.

Exudate containing these substances is physically bound within the dressing and safely removed during dressing change.





People.Health.Care.

Vliwasorb® sensitive border **NEW**

Superabsorbent wound dressing with silicone and adhesive border – the ideal addition to L&R's 3-step treatment recommendation.




Easy application
and atraumatic removal



Good exudate management



Reduction of germs

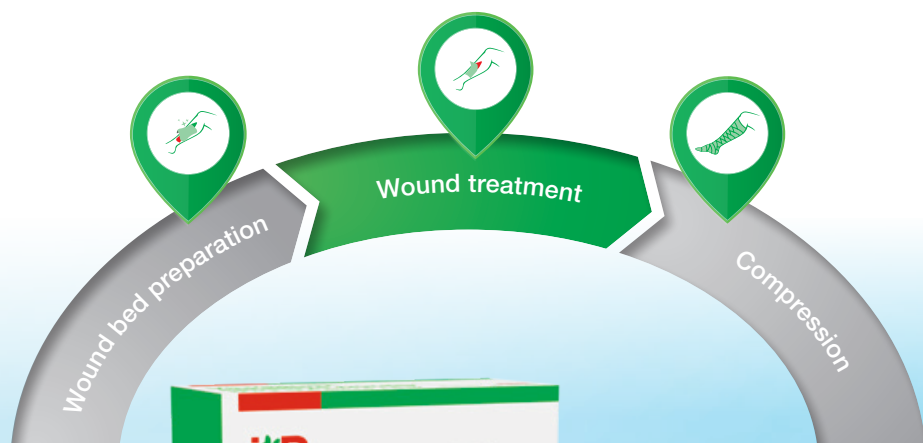

Ordering information

Vliwasorb® sensitive border **NEW**

sterile, individually sealed

Size Dressing (cm)	Size Wound pad (cm)	REF	Shipping units (pcs.)
10 x 10	5 x 5	159 558	10/200
12.5 x 12.5	7.5 x 7.5	159 559	10/100
17.5 x 17.5	12.5 x 12.5	159 560	10/100
15 x 25	9 x 19	159 561	10/100
20 x 25	14 x 19	159 562	10/100
Sacrum			
23 x 23	17 x 16.5	159 569	10/100

Publisher

Lohmann & Rauscher GmbH & Co. KG, Irlicher Straße 55, D-56567 Neuwied
Lohmann & Rauscher GmbH, Kirchengasse 17, A-2525 Schönaun an der Triesting

Contact addresses

Lohmann & Rauscher GmbH & Co. KG, Westerwaldstraße 4, D-56579 Rengsdorf
Lohmann & Rauscher GmbH, Johann-Schorsch-Gasse 4, A-1141 Wien
Lohmann & Rauscher AG, Oberstrasse 222, CH-9014 St. Gallen